

# PIN Authorization

3.24.2.0 | User Guide



**intelerad**

[Intelerad.com](https://www.intelerad.com)

## COPYRIGHT

© 2003-2024 DICOM Grid, Inc. dba Ambra Health. All Rights Reserved.

No portion of the contents of this publication may be reproduced or transmitted in any form or by any means without the express written permission of DICOM Grid, Inc. dba Ambra Health.

Information in this documentation is subject to change without notice and does not represent a commitment on the part of DICOM Grid, Inc. dba Ambra Health. The software described in this document is furnished under a license agreement or nondisclosure agreement and can be used and copied only in accordance with the terms of the agreement. No part of this document may be reproduced, or transmitted in any form or by any means, electronic or mechanical, including photocopying or recording, for any purpose, without the express written permission of DICOM Grid, Inc. dba Ambra Health.

Unless otherwise noted, names and clinical data used herein are fictitious. No association with any real person or clinical event is intended or should be inferred.

Except as expressly provided otherwise in writing, the information provided in this document is provided AS IS, without any condition or warranty, whether written, oral, express, implied, legal, or statutory. No condition or warranty of merchantability or fitness for a particular purpose applies to anything provided by Intelrad. Without limiting the foregoing, Intelrad nor any of its suppliers warrants this documentation is or will be accurate, error free or meets or will meet user requirements.

## TRADEMARKS

Intelrad®, Intelrad Medical Systems, and the Intelrad logo, IntelePACS, InteleViewer™, Enterprise Viewer™, InteleOne®, InteleArchive™, InteleShare VNA™, InteleShare™, InteleShare Research™, InteleOrchestrator™, InteleHeart™, InteleScreen™, and InSight PACS™ are either registered trademarks or trademarks of Intelrad Medical Systems Incorporated or its affiliates.

## THIRD-PARTY TRADEMARKS

Adobe, Acrobat, and Reader are either registered trademarks or trademarks of Adobe Systems Incorporated in the United States and/or other countries. Firefox is a registered trademark of Mozilla Foundation in the United States and other countries. Google Chrome browser is a registered trademark of Google Inc. in the United States and other countries. Internet Explorer is a registered trademark of Microsoft Corporation in the United States and other countries. Intel, Pentium, Pentium II Xeon, and Pentium III Xeon are trademarks or registered trademarks of Intel Corporation or its subsidiaries in the United States and other countries. Microsoft and Windows are either registered trademarks or trademarks of Microsoft Corporation in the United States and/or other countries. NVIDIA is a trademark or registered trademark of NVIDIA Corporation in the United States and/or other countries. Oracle and Java are registered trademarks of Oracle Corporation and or its affiliates. PowerScribe 360 is a trademark or registered trademark of Nuance Communications Inc. or its affiliates in the United States and/or other countries. Sun, Sun Microsystems, and Java are trademarks or registered trademarks of Sun Microsystems, Inc. in the United States and other countries. Safari, Mac, and OS X are trademarks of Apple Inc., registered in the United States and other countries.

All other brand names, product names, or trademarks belong to their respective holders.

## INTENDED USE

InteleShare PACS software is intended for use as a primary diagnostic and analysis tool for diagnostic images for hospitals, imaging centers, radiologists, reading practices and any user who requires and is granted access to patient image, demographic and report information.

InteleShare ProViewer, a component of InteleShare PACS, displays, modifies and manages diagnostic quality DICOM images including 3D visualization and reordering functionality.

Lossy compressed mammographic images and digitized film screen images must not be reviewed for primary image diagnosis or image interpretations. Mammographic images may only be viewed using cleared monitors intended for mammography display.

Not intended for diagnostic use on mobile devices.

CONTRAINDICATIONS—None.

Caution: Federal law restricts this device to sale by or on the order of a physician. Fluency in English is required as a prerequisite for operating the Ambra PACS product.

This system does not replace the education, skill, and judgment of properly trained medical practitioners. Only properly trained and qualified individuals shall have access to and use InteleShare PACS and must know of its functionality, capabilities and limitations. Typical users of this system are trained health professionals, physicians, nurses, and technologists.

Downloaded Images, Workstations and Isolated Installs: You and your users must maintain InteleShare PACS with the most current versions, including available updates and upgrades. Delaying or refusing updates or upgrades following a recall may result in a non-compliant InteleSharePACS.

SAFETY ISSUES: InteleSharePACS is a medical device, and as such, must meet medical device safety and effectiveness requirements imposed by national regulations. Any unmonitored or unconnected use of InteleSharePACS, or use of InteleSharePACS without a valid right may put the health and safety of patients at risk as you will not be advised of the availability of any software patch, bug fix, update or upgrade nor will be informed of Field Safety Notices, Medical Device Recalls or Advisory Notices related to InteleSharePACS. Client and authorized users must consult national regulatory site(s) to be informed of Field Safety Notices, Medical Device Recalls or Advisory Notices related to InteleSharePACS. Intelrad does not have access to authorized users systems to implement corrections to prevent (or correct) occurrences of patient safety issues. You are responsible to flow down recall and patient safety information to your users. The user of the medical device must report any serious incident that has occurred in relation to the medical device to the manufacturer (Intelrad) and the competent authority having jurisdiction in their locale.

Referring Physicians Use: Images for authorized referring physicians may not be of diagnosis quality and should not be used for diagnostic purposes.



DICOM Grid, Inc. dba Ambra Health  
305 Church at N Hills St, 6th Floor,  
Raleigh, North Carolina 27609, USA



InteleShare UDI is: **+AMBRAPACS0/\$73.24.2.0D**

## DECLARATION OF CONFORMITY

Declaration of Conformity is issued under the sole responsibility of the DICOM Grid, Inc. dba Ambra Health.

Title: PIN Authorization User Guide

Software version: 3.24.2.0 (Multiple Patient Multiple Use)

Date: 2024-05-17

Part number: AMPINEN3.24.2.0UG-O Issue 001

# TABLE OF CONTENTS

<b>Browser Requirements</b> .....	<b>4</b>
<b>Contacting Intelrad Technical Support</b> .....	<b>5</b>
<b>Obtaining Printed Documentation</b> .....	<b>7</b>
<b>PIN Authorization</b> .....	<b>8</b>

# 3

## BROWSER REQUIREMENTS

InteleShare supports the following browsers:

Microsoft Internet Explorer 9 or later

Microsoft Edge

Apple Safari

Google Chrome

Firefox.

# 4

## CONTACTING INTELERAD TECHNICAL SUPPORT

Your PACS administrator can assist you with any issues you may encounter. If you require additional assistance, you can contact IntelShare Technical Support, 24 hours a day, seven days a week.

To contact us:	Use:
On the Internet	<a href="https://serviceportal.intelerad.com/csm">https://serviceportal.intelerad.com/csm</a>
By telephone	Toll-free: 888-315-0790

These coordinates and a wealth of other information are also available on the Intelerad Service Portal.

<https://serviceportal.intelerad.com/csm>

You should regularly check the Intelerad knowledge base for the latest version of the documentation.

IntelShare PACS is a cloud-based product. Customers are informed of upcoming releases via multiple channels to ensure they stay up-to-date with the latest developments and enhancements. All customers receive timely and comprehensive e-mails from our marketing team, which not only highlight the release dates but also describe in detail the new features and functionalities that will be available in the upcoming release. When a new version of the software becomes available it will automatically upgrade on the appointed release date of that version. Please contact support with for any questions regarding upgrades.

When you contact IntelShare Technical Support to report a problem, please have at hand the following information, as applicable:

- client code and location of your IntelShare PACS installation
- full error message and the steps required to reproduce the problem
- AE Titles of the affected devices
- operating systems of any affected machines
- description of the problem and when it first occurred

If the problem affects a particular study, please also provide the following:

- patient ID or patient number (M.R.N.)
- accession number/requisition number
- modality type and name

For security information, please open a ticket with Intelrad Technical Support.

# 5

## OBTAINING PRINTED DOCUMENTATION


Intelerad offers printed and bound versions of product documentation free of charge. To request printed copies of Intelerad documentation, contact your Client Success manager. The printed documents will be provided within 7 days or less.

# 6

## PIN AUTHORIZATION

Click the link in your EMR to access your study images. A 'Sign In' window will open.

Enter your email address and click 'View Study'.

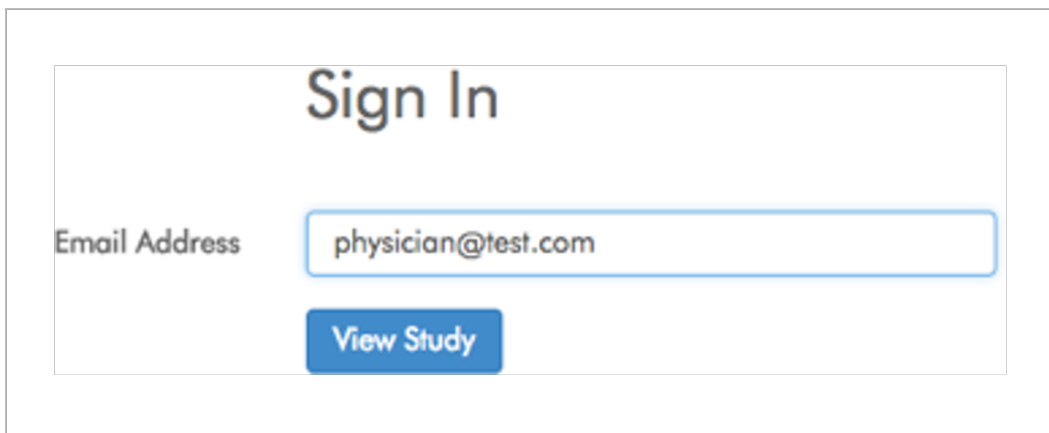


The image shows a 'Sign In' form. At the top, the text 'Sign In' is displayed in a large, dark font. Below this, there is a label 'Email Address' on the left. To the right of the label is a text input field containing the email address 'physician@test.com'. Below the input field is a blue button with the text 'View Study' in white.

A PIN number will be sent to your email address.

Enter the PIN number and click 'Submit'.

Your study images will appear.



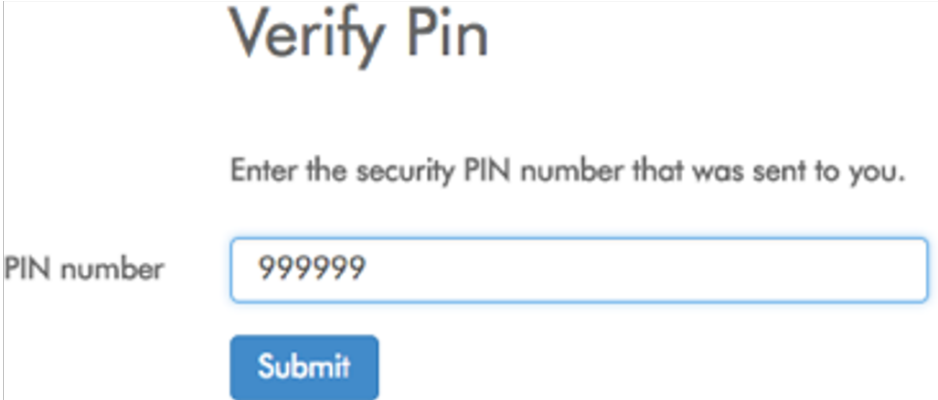

The image shows a 'Verify Pin' form. At the top, the title 'Verify Pin' is displayed in a large, bold, dark font. Below the title, a subtitle reads 'Enter the security PIN number that was sent to you.' To the left of the input field, the label 'PIN number' is shown. The input field itself is a rounded rectangle with a blue border, containing the text '999999'. Below the input field is a blue rectangular button with the word 'Submit' written in white text.

- 📋 Your email address must be registered with the organization. If it is not, you will be alerted to contact your administrator.
- 📋 PINs expire after 10 minutes.
- 📋 A PIN entered incorrectly three times is invalidated.